



A. Setup


1 Humidification chamber preparation

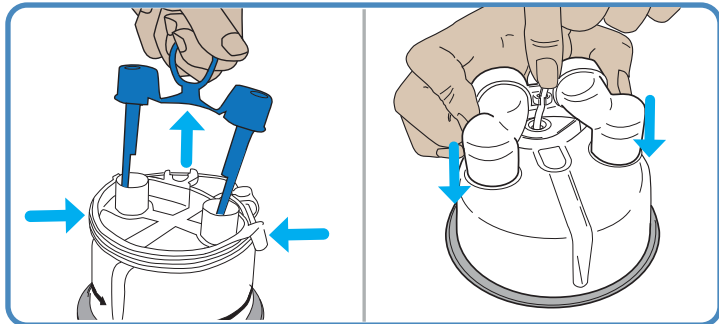


2 Chamber fitting

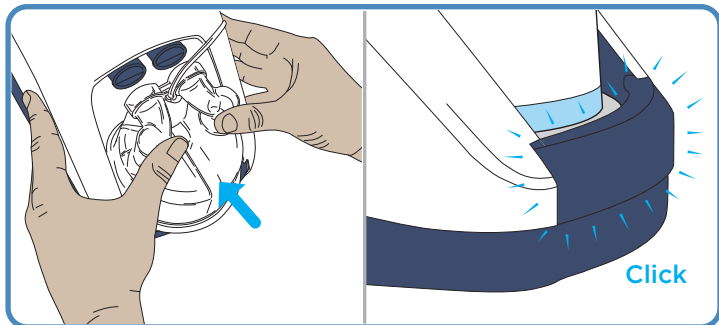

3 Water bag connection

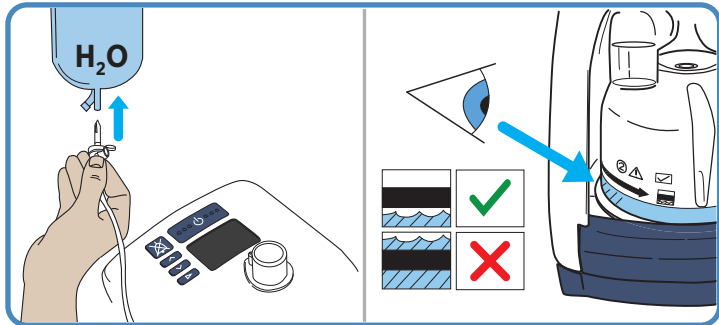

4 Breathing tube connection

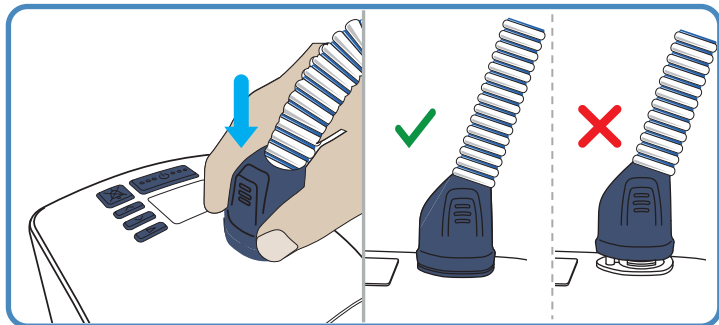

Scan the relevant QR code above to download the Airvo 2 simulator app to your smartphone.

Fisher & Paykel
HEALTHCARE

B. Start


1 Press 5 seconds




GREEN LIGHT:
Correct disinfection has
been completed since
last use



ORANGE LIGHT:
Disinfection
NOT completed

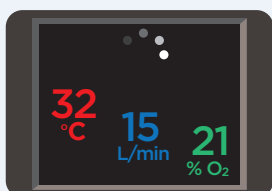


ROLLING CIRCLE:
Warm up time

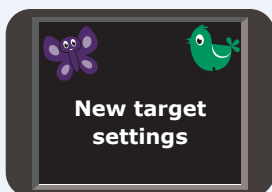

2 To activate/deactivate JUNIOR MODE:

Press 5 seconds

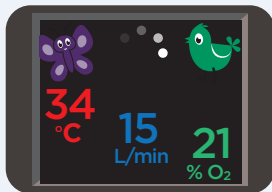


3 For Target Flow and % O₂ settings see **C**


C. Settings

1 Target Flow Setting

Access the target flow setting.

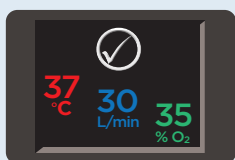


Move past the temperature setting menu.

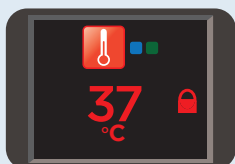


Adjust target flow:



If the setting menu is locked (🔒), unlock by pressing 3 sec: ⏏

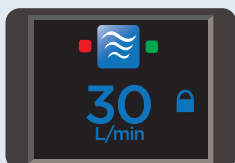

Parameter setting adjustment

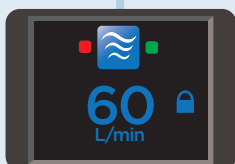


Exit the menu settings.




2 % O₂ Setting

Adjust the O₂ flow meter to change % O₂ value.

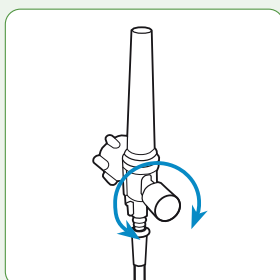

Note: This will NOT affect the delivered patient flow.

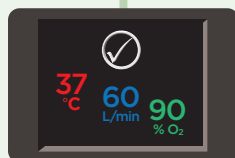

Read the measured % O₂ value.


3 When the target flow setting is adjusted:

New **target flow setting** adjustment (item 1) will change the **delivered % O₂**.


If required, adjust the % O₂ (item 2).

Patient Interfaces

Nasal Interface

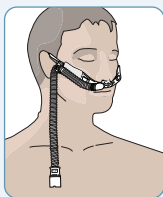


OPT94X
OPT96X

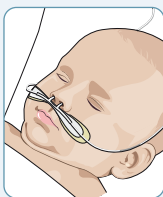


OJR416
OJR418
OJR520

Mask Interface

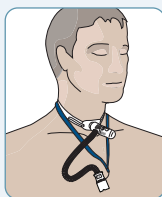

OPT970

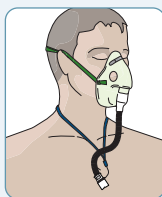


OPT980

Available Target Settings

		Temp °C			Flow L/min						
		31	34	37	2	5	10	15	20	25 - 55	60
OJR416	L		●		2 20						
OJR418	XL		●		2 25						
OJR520			●	●	10 50						
OPT942/ OPT962	S		●	●	10 50						
OPT944/ OPT964	M		●	●	10 60						
OPT946/ OPT966	L		●	●	10 60						
OPT970				●	10 60						
OPT980		●	●	●	10 60						


Turn off

Press 5 seconds

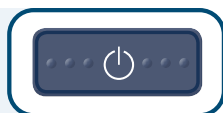


Cleaning and Disinfection

A. Clean the outlet elbow

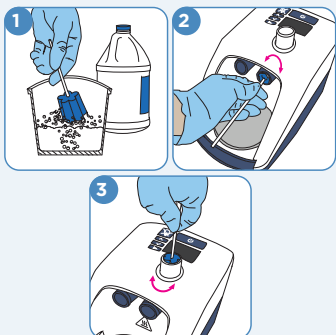

B. Surface cleaning

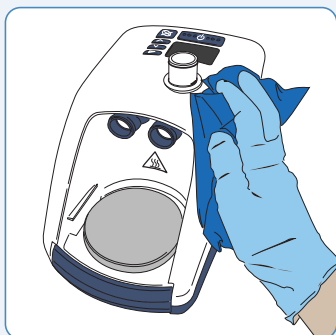

C. Automatic disinfection cycle

